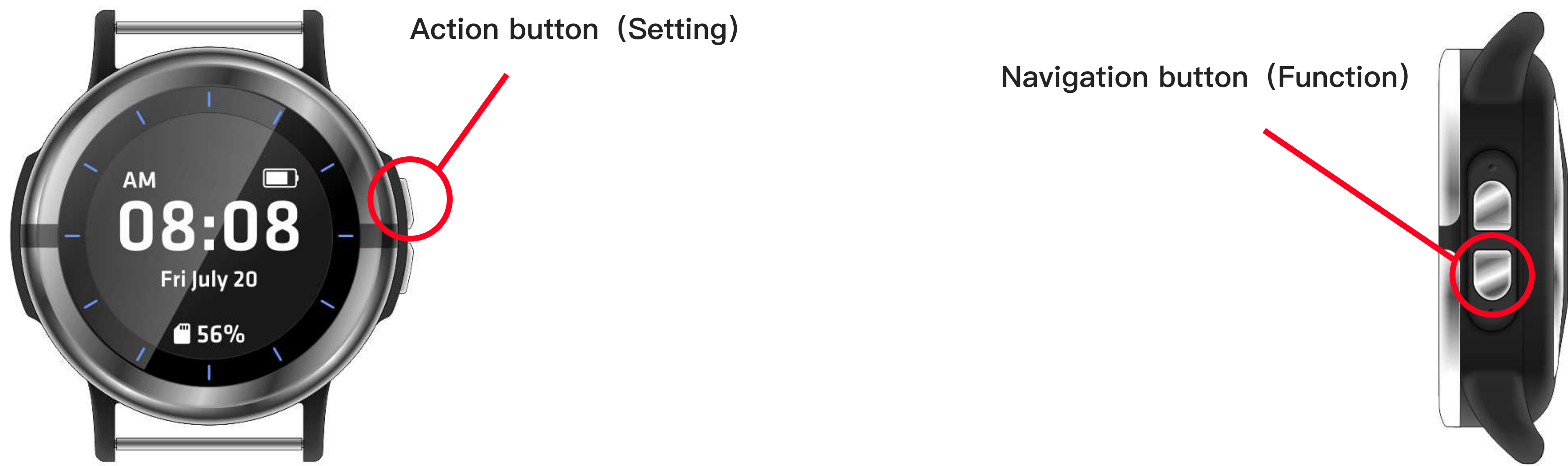


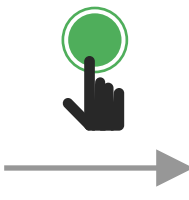
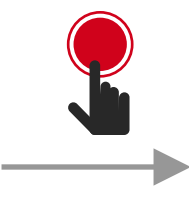
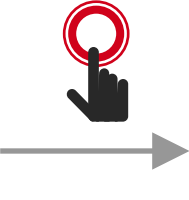
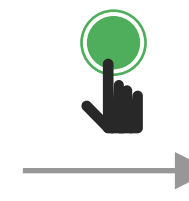

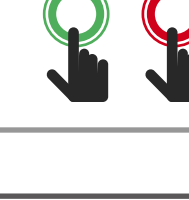




UX Version: ADI\_UX\_210115\_V2.2


Designer :Holiday




Gestures specifications

Function page	Navigation button		Navigation function
	Action button		Confirm
	Action button Long press		Quit/Back
In shutdown mode	Navigation button		Starting up
work mode	Navigation button Long press		Restart
work mode	Navigation button + Action button Long press		Upgrade
Shipment mode	USB		Charging activation
			Automatic jump

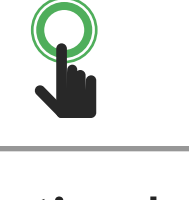
Gesture Specifications for Automatic shutdown

Automatic shutdown	BLE Disconnected	static page	dynamic page	power on
	USB cable Disconnected	30s+1min	60s+1min	 Navigation button Long press
	Filesystem Logging stopped			
	No key press detected			

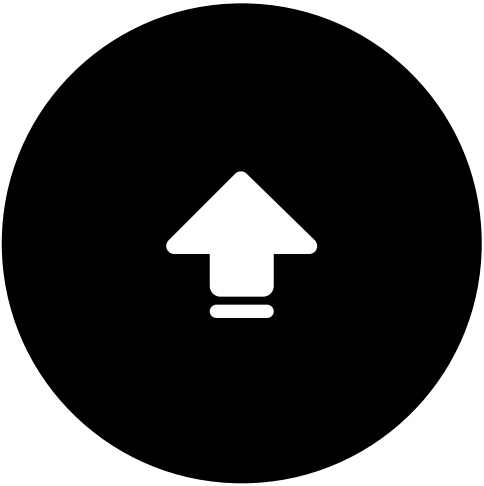
Dynamic Page




updating

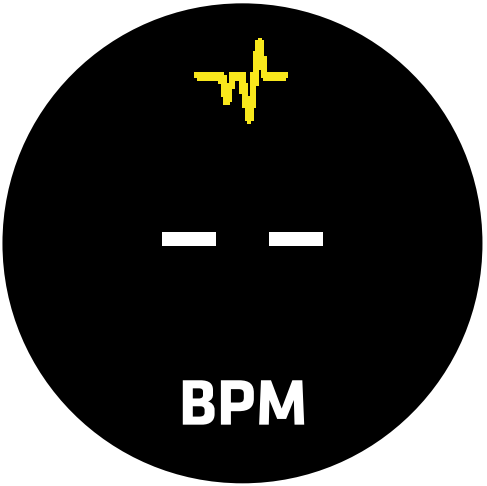


BPM

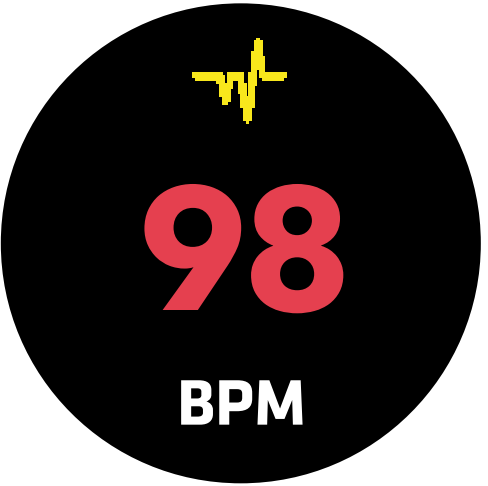


98  
BPM

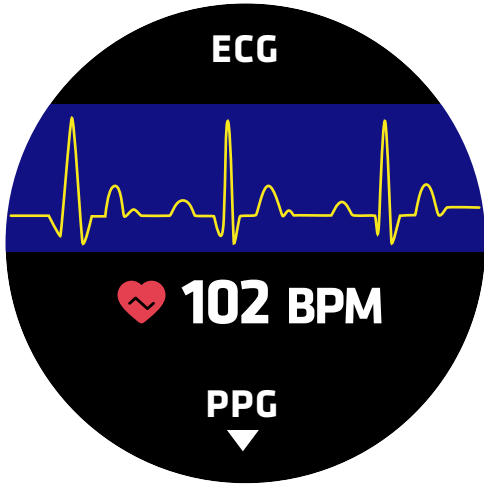

ECG  
102 BPM  
PPG

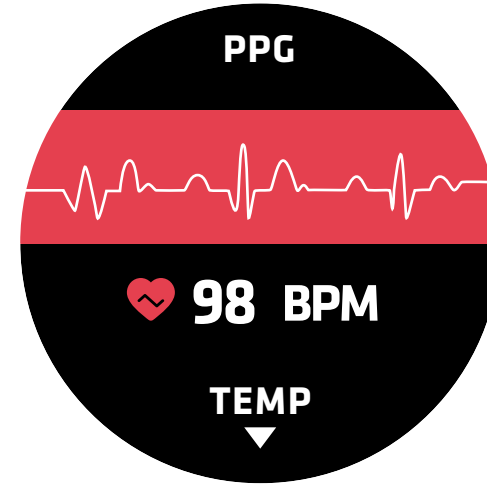


PPG  
98 BPM  
TEMP

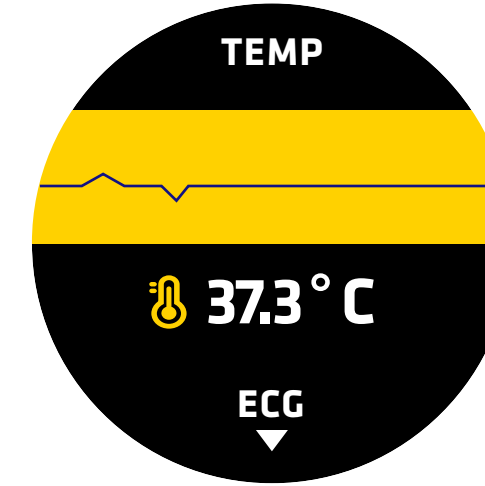

TEMP  
37.3 °C  
ECG


Starting up、Restart、Uparade

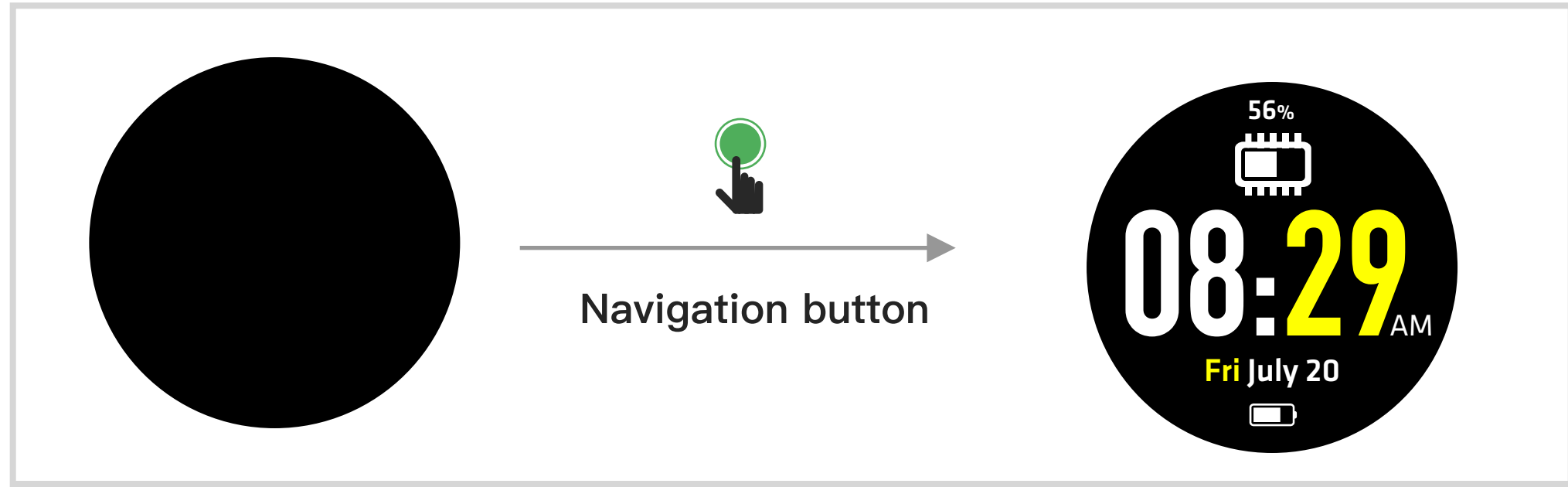starting up




In power off mode, press the Navigation button for 1s to start the machine

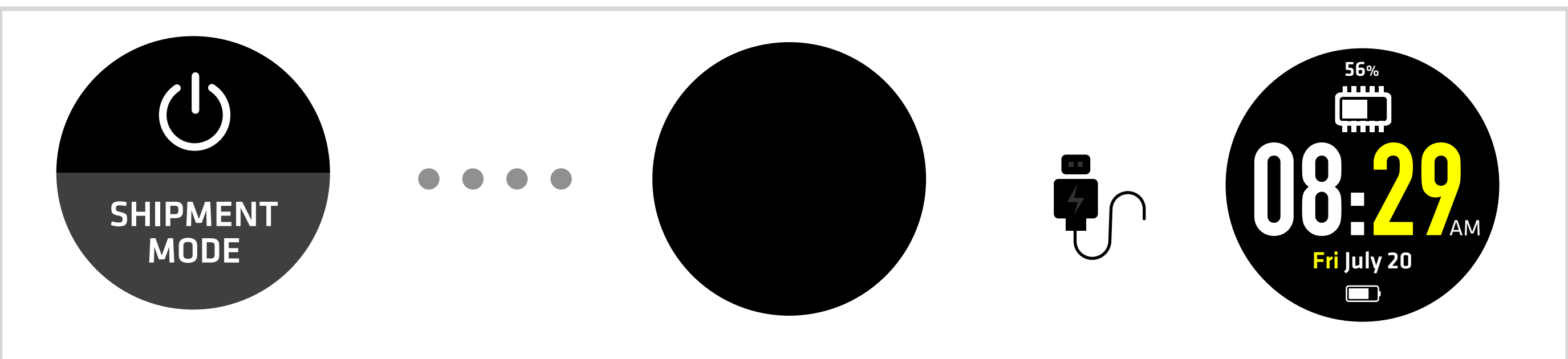starting up in shipment mode




It must be started with a charging cable

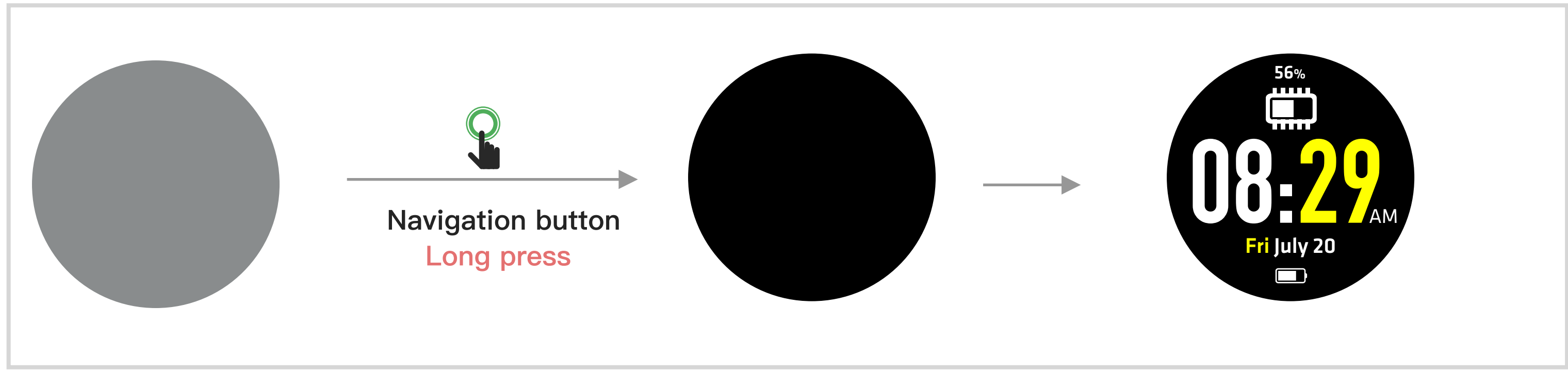
restart



upgrade

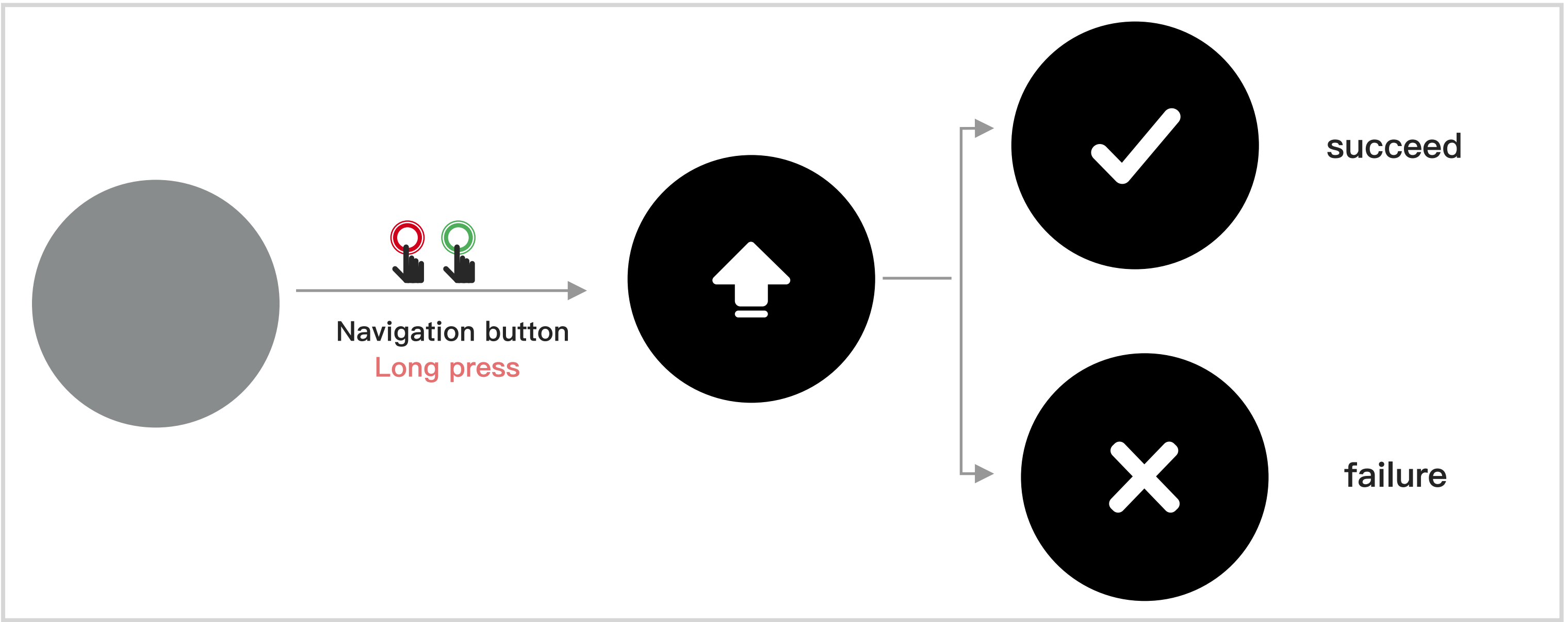


Main interface (static page)

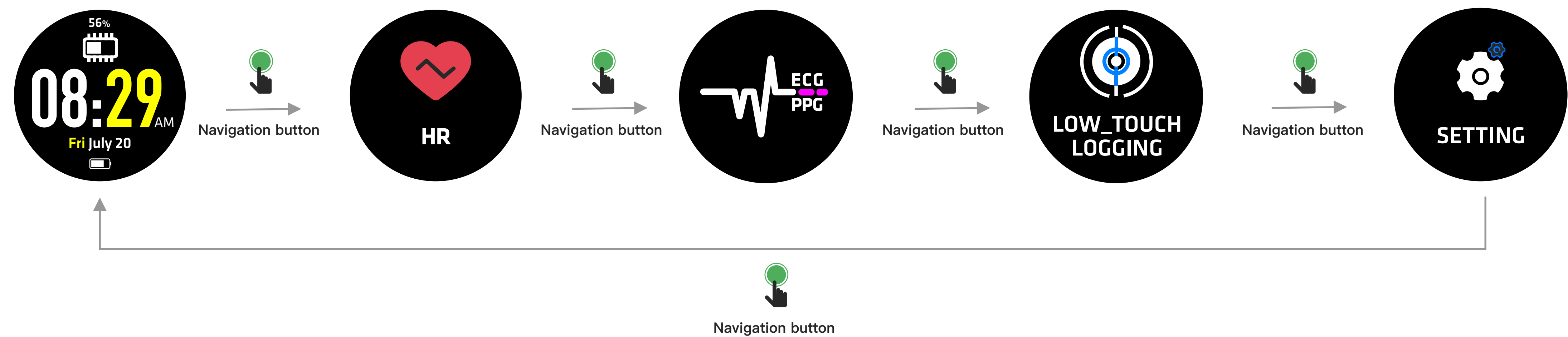

 short press: Navigation function

No action in current Page: Backlight will auto close after 30s

# Heart rate

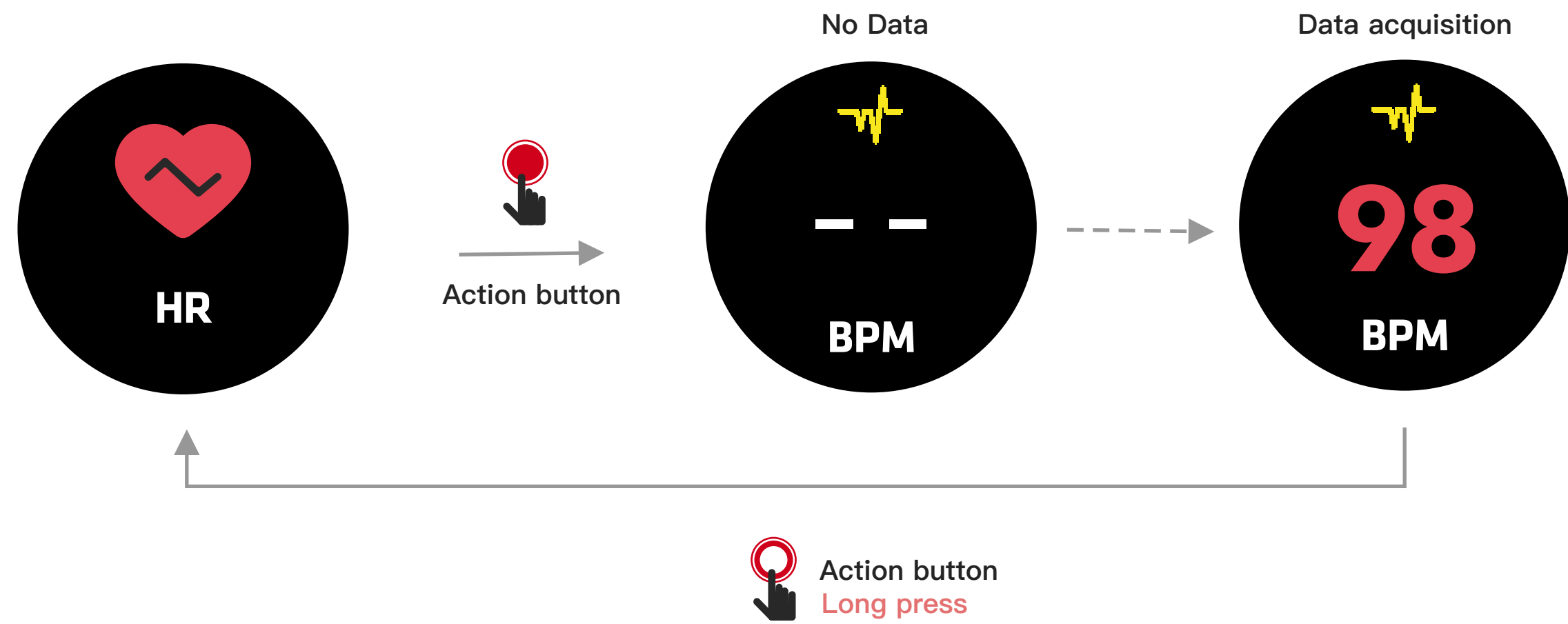


short press: Confirm




long press: Quit

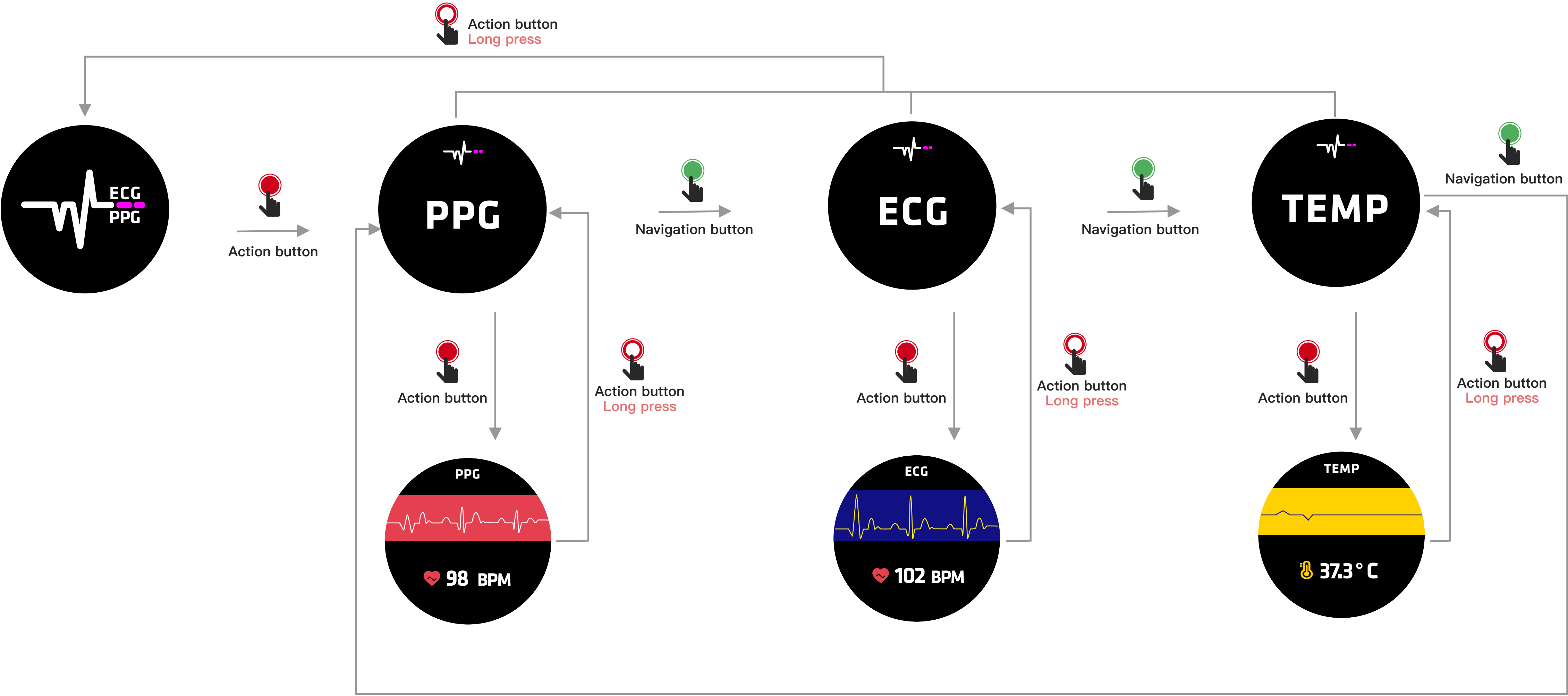Waveform




No action in current Page: Backlight will auto close after 30s

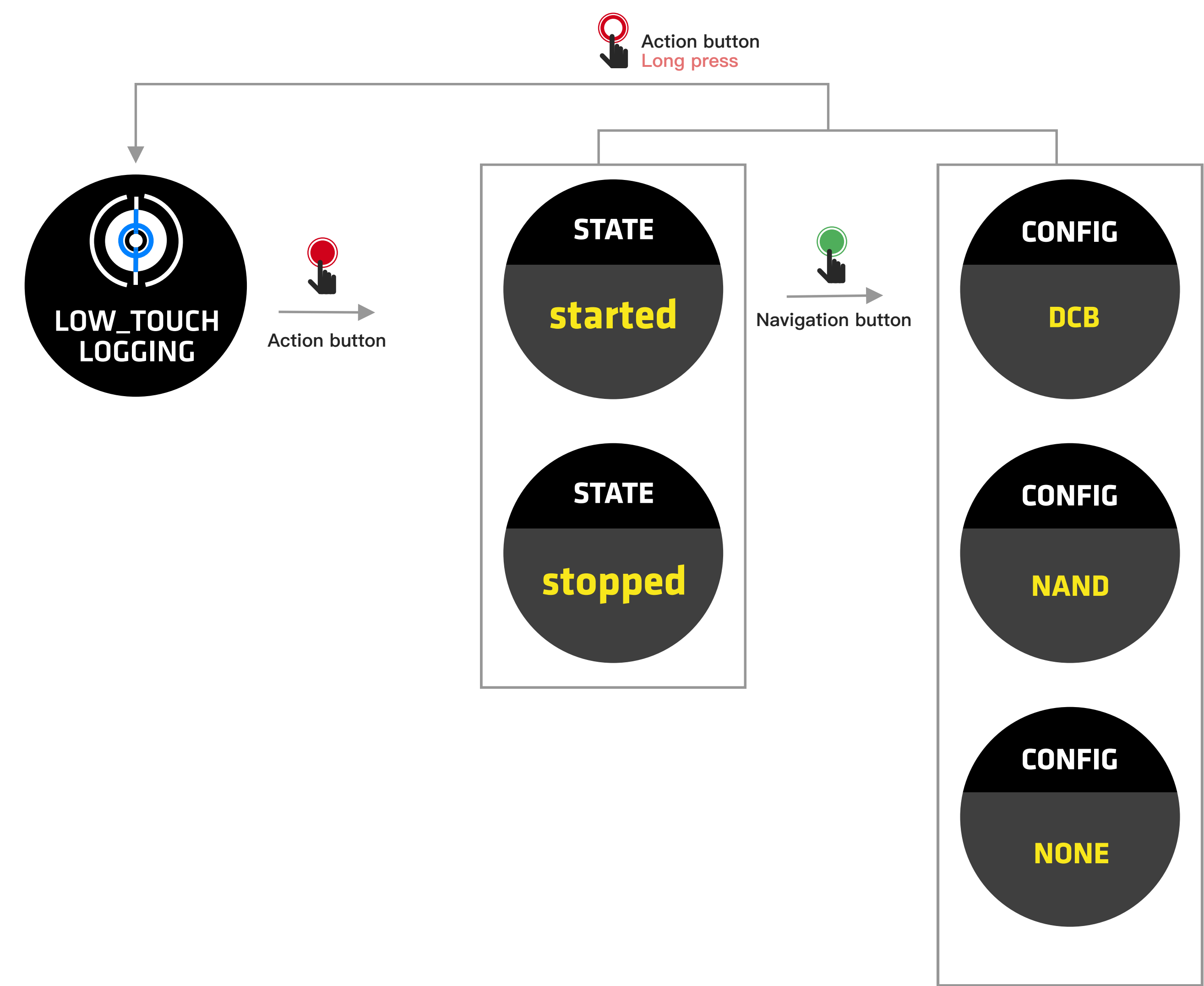LOW\_TOUCH LOGGING



 short press: Confirm


 long press: Quit/Back

 short press: Navigation function

No action in current Page: Backlight will auto close after 30s




## Setting (static page)

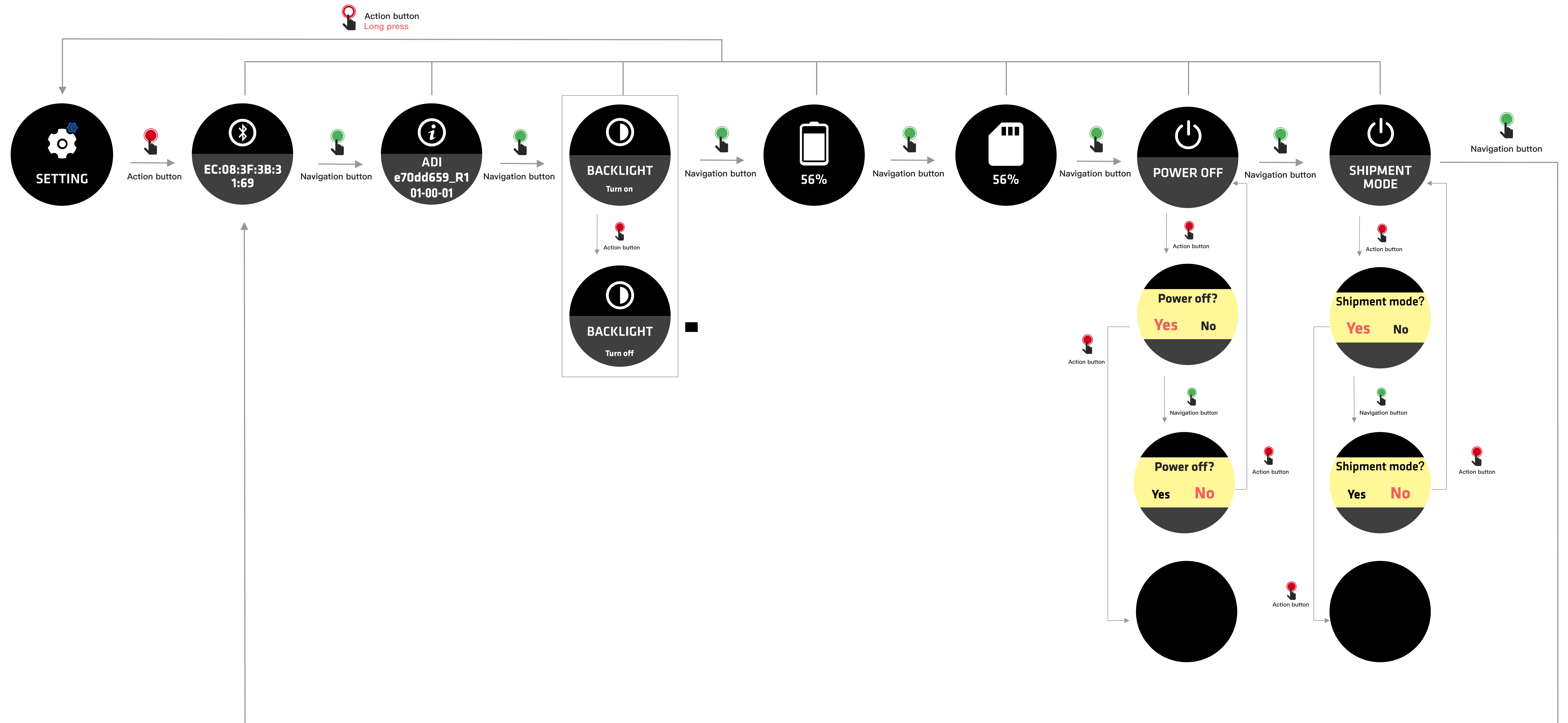

 short press: Confirm

 long press: Quit/Back

 short press: Navigation function

**No action in current Page: Backlight will auto close after 30s**